

Oinsá bele halo Teste Antígenu Rápidu (RAT)

KEEP IT
COVID
safe

Iha estoju RAT oioin ne'ebé disponivel.
Tenke halo tuir instrusaun ne'ebé fó ona iha estoju teste ne'e.

Tenke halo tuir matáðalan tuir mai:

1. Lee didi'ak instrusaun iha estoju laran molok ita hahú.
Instrusaun balu iha teste inklui kódiçu QR ne'ebé Ita bele scan hodi haree vídeo.




2. Fase liman.



3. Loke embalajen.



4. Dulas swab iha inus laran ho kle'an 2cm ho movimentu sirkulár. Dala lima ka segundu 15 iha kada inus kuak.

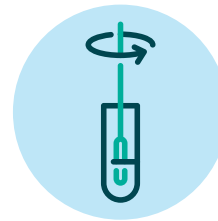

5. Tau swab iha tubu solusaun. Kahur didi'ak. Dala ruma tenke ku'u swab ka fera stick - haree instrusaun iha teste.




6. Taka tubu ho ninia matan.



7. Tau bee-musan iha janela amostra husi teste. Ita-nia instrusaun iha estoju teste sei hatete bee-musan hira. Labele book ka hiit estoju teste.



8. Hahú relójiu-programadór. Heín durante tempu ne'ebé hatuur ona iha instrusaun. Labele book ka hiit Janela Amostra.

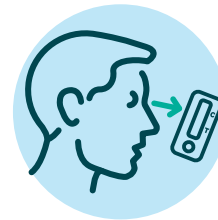

9. Haree rezultadu (haree ezemplu iha liman-loos).

Reading your result



'C' Kontrolu (líña iha ne'e hatudu katak teste funsiona)
'T' Teste (Iaiha líña = negativu, líña iha ne'e = pozitivu)




NEGATIVU

Iha possibilidade ki'ik katak ita hetan COVID.
Sekarik Ita iha sintoma, repete teste RAT ka halo teste PCR iha loron tuir mai.



INVALIDU

Halo tan teste (no repete pasu 2-9). Iha difikuldade? Halo teste PCR.

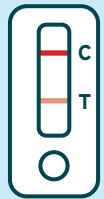

POZITIVU

Ita hetan COVID. Karik rezultadu teste iha 'T' hatudu líña, mezmuke mihis, ida-ne'e rezultadu pozitivu.

Se Ita hetan rezultadu pozitivu, izola an kedas no fó-hatene Public Health.

Scan kódiçu ne'e hodi notifika rezultadu online ka telefona Public Health Hotline 1800 671 738.

